



Aurora Masks by Inogen®


Aurora F1 Full Face Mask

The Aurora F1 full face mask delivers comfort and performance—designed to support ease of adoption and long-term adherence.

In a 30-day paired analysis (N=46 full face mask), patients reported higher satisfaction with the Aurora F1 compared to their current mask.¹



Aurora F1 Full Face Mask

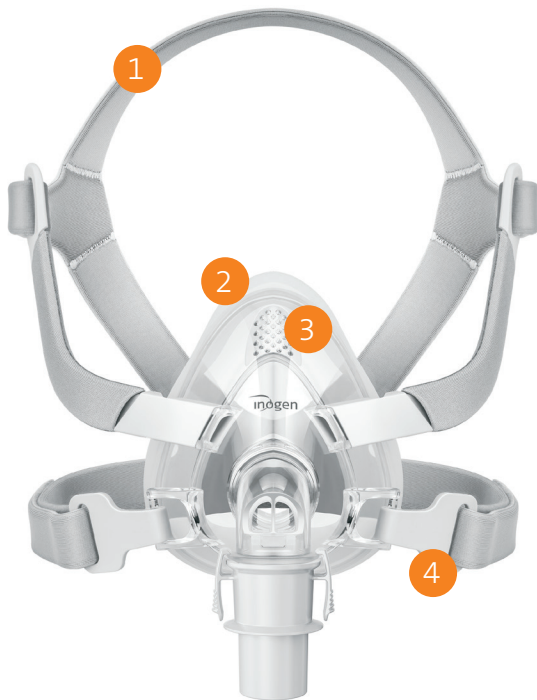

SPECIFICATIONS

Therapy Pressure Range	4-30 cmH ₂ O
Tube Connection	22mm
Flow Resistance (At 50 L/min)	≤0.5 cmH ₂ O
Flow Resistance (At 100 L/min)	≤1.5 cmH ₂ O
Expiratory Resistance with Anti Asphyxia Valve Open (At 50 L/min)	≤2.0 cmH ₂ O
Inspiratory Resistance with Anti Asphyxia Valve Open (At 50 L/min)	≤2.0 cmH ₂ O
Anti Asphyxia Valve Open Pressure	≤2.0 cmH ₂ O
Anti Asphyxia Valve Closed Pressure	≤2.0 cmH ₂ O
Cushion Sizes	S/M/L/XL

In a 30-day study (N=46 full face mask users), 94% found the Aurora F1 mask to be just as or more durable as their current mask.¹

1 High-Quality Ergonomic Crown-Shape Headgear

Mask headgear design features a triple-material design with ultrasonic edge sealing technology, promoting comfort and durability.

2 Adaptive Silicone Cushion

Medical-grade silicone cushion design features an adaptive design with a range of up to 15 degrees at the nose bridge, adjusting to different face shapes.

3 Quiet Honeycomb Vent Design

Honeycomb vent design reduces airflow intensity to provide gentle and quiet exhalation for patients, designed to improve comfort and adherence.

4 Quick-Release Headgear Clips

Quick-release headgear clips allow for ease of putting on and taking off the mask, making adjustments easy and user-friendly.

provider.inogen.com

1. AMG CPAP Mask Trial Study Aug 2025, 30-day study of 46 patients using CPAP therapy with a full face mask for: at least 6 months, at least 5 nights a week, and at least 4 hours a night.

Inogen® Aurora Masks

The mask is a noninvasive accessory used for channeling airflow (With or without supplemental oxygen) to a patient from a positive airway pressure (PAP) device such as a continuous positive airway pressure (CPAP) or bi-level system. For safe and proper use of the Inogen® Aurora masks (including intended use, contraindications, and warnings), please refer to the instructions for use at www.inogen.com/auroramasks.

PM250182 EN_EX_USA | Rx Only. © 2025 Inogen, Inc. | 859 Ward Drive, Suite 200, Goleta, CA 93111

Inogen® and other related product names contained herein are trademarks of Inogen, Inc. The usage of any Inogen, Inc. trademark is strictly forbidden without the prior consent of Inogen, Inc. All other trademarks are trademarks of their respective owners or holders.